Instructions for Use



Important

Please read the safety information and all information delivered with the product carefully to familiarize yourself with safe and effective usage.



Legal information

Warning notice system

This manual contains notices you have to observe in order to ensure your personal safety, as well as to prevent damage to property. The notices referring to your personal safety are highlighted in the manual by a safety alert symbol, notices referring only to property damage have no safety alert symbol. These notices shown below are graded according to the degree of danger.

↑ DANGER

indicates that death or severe personal injury will result if proper precautions are not taken.

⚠ WARNING

indicates that death or severe personal injury may result if proper precautions are not taken.

↑ CAUTION

indicates that minor personal injury can result if proper precautions are not taken.

NOTICE

indicates that material damage can result if proper precautions are not taken.

If more than one degree of danger is present, the warning notice representing the highest degree of danger will be used. A notice warning of injury to persons with a safety alert symbol may also include a warning relating to property damage.

Qualified personnel

The product/system described in this documentation may be operated only by **personnel qualified** for the specific task in accordance with the relevant documentation, in particular its warning notices and safety instructions. Qualified personnel are those who, based on their training and experience, are capable of identifying risks and avoiding potential hazards when working with these products/systems.

Use of EIZO products

↑ WARNING

EIZO products may only be used for the applications described in the catalog and in the relevant technical documentation. If products and components from other manufacturers are used, these must be recommended or approved by EIZO. Proper transport, storage, installation, assembly, commissioning, operation and maintenance are required to ensure that the products operate safely and without any problems. The permissible ambient conditions must be complied with. The information in the relevant documentation must be observed.

Trademarks

All names identified by ® are registered trademarks of their respective owners. Please refer to the trademarks listed in the appendix. The remaining trademarks in this publication may be trademarks whose use by third parties for their own purposes could violate the rights of the owner.

Disclaimer of liability

We have reviewed the contents of this publication to ensure consistency with the hardware and software described. Since variance cannot be precluded entirely, we cannot guarantee full consistency. However, the information in this publication is reviewed regularly and any necessary corrections are included in subsequent editions.

Table of contents

	Legal information						
1	Intro	ductionduction	. 5				
	1.1	Contents of this document	5				
	1.2	Intended use	5				
	1.3	User	6				
2 Safety information							
	2.1	General safety instructions	7				
	2.2	Product-specific safety information	11				
3	Desc	eription	12				
	3.1	Scope of delivery	12				
	3.2	Monitor performance features	12				
4	Setup	p and installation	14				
	4.1	Installation site	14				
	4.2	Installing the monitor	16				
5	Conn	necting	17				
	5.1	Safety information for connection	17				
	5.2	Connector locations					
	5.3	Connection panel	19				
	5.4	Connecting the signal cable	21				
	5.5	Connecting the power cable	22				
6 Commissioning							
	6.1	Switching on the monitor and video source					
	6.2	Avoiding image sticking					
	6.3	Check for pixel defects					
	6.4	Setting the image geometry					
	6.5	Adjustment of monitor and video source					
	6.6	Adapting the monitor using Force Mode	26				
		6.6.1 Introduction to Force Mode	26				
		6.6.2 Basic information on timing					
		6.6.3 Entering timing data					
_		6.6.4 Determining the timing data					
7	•	ration					
	7.1	Operator controls					
	7.2 Remote control						
	7.3	Lock or unlock OSD menu					
	7.4	Picture layout (PaP, PiP, PoP)					
	7.5	Overview of OSD menu	37				

	7.6	"Picture" menu	37
	7.7	"Image" menu	39
	7.8	"Signal" menu	42
	7.9	"LUT" menu	46
	7.10	"Info" menu	47
	7.11	"Utilities" menu	48
	7.12	"Force Mode" menu	50
	7.13	"Adjust LUT" menu	51
8	Clean	ing and check settings	53
	8.1	Cleaning	53
	8.2	Check the settings	54
9	Troub	eleshooting	55
	9.1	Fault correction	55
10	Techr	nical specifications	56
	10.1	Monitor characteristics	56
	10.2	Power supply	56
	10.3	Electronics	57
	10.4	Mechanical design	57
	10.5	Climatic conditions	58
	10.6	Safety regulations	58
11	Dime	nsion drawings	59
	11.1	View from front and side	59
	11.2	Rear view	59
12	Spare	parts / accessories	60
	12.1	Accessories	60
13	Appe	ndix	61
	13.1	Information on electromagnetic compatibility (EMC)	
	13.2	Markings and symbols	66
	13.3	FCC Declaration of Conformity	67
	13.4	China RoHS (Restriction of Hazardous Substances)	68
	13.5	Declaration of compliance with India RoHS	69
	13.6	Environmental protection	69
	13.7	Additional devices	69
	13.8	Warranty	69
	13.9	Repairs	70
	13.10	Contact	70
	13.11	Trademarks	70
	Index		71

1 Introduction

1.1 Contents of this document

This document explains the functionality and the approved use of the CuratOR LX491W. To ensure clarity, it does not contain all detailed information on this product.

The contents of this document are neither part of a previous or existing agreement, commitment or legal relationship, nor does it modify such.

Note

This documentation is available in electronic format only. It is included on the CD-ROM provided and can be downloaded at www.eizo-or.com, or provided by the sales partner from whom you purchased the product.

1.2 Intended use

Intended purpose

The CuratOR LX491W is intended for the display of still images and moving images from various commercially available devices commonly used in a medical environment, in particular endoscopic. The monitor is optimized for the reproduction of colour images. The monitor is not suitable for mammography.

Intended patient population and medical conditions

The LX491W can be used for the intended purpose irrespective of age, body weight and gender.

The LX491W is intended to be used in combination with medical devices. The monitor therefore has no direct contact with the patient.

The LX491W is intended to display still images and moving images from various commercially available (medical) devices commonly used in a medical environment. The monitor cannot be used for direct diagnosis and as main device for monitoring live support equipment.

Intended users

The intended users for the LX491W are qualified healthcare professionals.

Intended environment

The LX491W is intended to be used in professional healthcare facilities such as clinics and hospitals. The monitor can be used in operating rooms (OR) or near patients, but is not limited to them. The monitor is not intended for direct patient contact!

The LX491W is not suited for the following environments:

- Home-based healthcare facilities.
- Near short-wave therapy devices.
- Near an MRI-System.
- · Built into vehicles, including ambulances.

1.3 User

Note

Serious incident

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

1.3 User

User

In the following, healthcare personnel such as surgeons or medical technicians are referred to as the "user".

Service / service personnel

"Service" or "Service personnel" identifies authorized personnel with knowledge of electrical and signal connection, local standards for image quality requirements, and safety of medical products, for example a hospital technician or manufacturer of medical devices.

Cleaning staff

"Cleaning staff" refers to personnel responsible for cleaning medical devices.

2 Safety information

2.1 General safety instructions

Correct and safe operation of EIZO devices assume professional transport, storage, installation, and connection, as well as careful operation and service.

The devices may only be used for applications for which they are intended.

For safety reasons, the following precautions must be observed:

DANGER

Please observe all warning information present on the device and in the instructions for use.

There is a danger to life if warnings are not obeyed. Severe personal injury or damage to property may occur.

Observe the safety requirements of EN 60601-1 (IEC 60601-1)

To prevent injury to patients and users, connect the electrical system in accordance with the safety requirements of EN 60601-1 (IEC 60601-1) for "Safety requirements for medical electrical systems".

Connecting the protective earth conductor

If the device is connected to line power, the device must be connected to a protective ground conductor. This is the only way to ensure that the touch leakage current in a first fault event does not exceed 500 μ A.

The interruption of the device's protective conductor is considered a first fault event in accordance with EN 60601-1.

Use the following measures to ensure that the leakage currents remain below the specified limits:

- Separators for signal input unit or signal output unit
- · Use of a safety isolating transformer
- · Use of the additional protective ground terminal

Mounting of the monitor: The monitor's suspension arm must have its own protective ground conductor. This protective ground conductor guarantees, together with the protective ground conductor of the monitor, that the housing leakage current always remains less than 500 μ A, even in the event of a single fault condition.

No unauthorized opening of the device / no unauthorized service or maintenance work

The device may only be opened by qualified personnel. Likewise, service or maintenance work may only be carried out by qualified personnel. There is a risk of electric shock.

No liability is accepted for death and injury to persons or damage to property resulting from work carried out by non-qualified personnel.

Do not touch components in the device

If the device is connected to the line power, components in the device are subjected to high voltages. Touching the components may be fatal.

No contact between device and patients

The device is not suitable for direct contact with a patient. The device and patient must never be touched simultaneously. Otherwise there is a danger to life and limb.

2.1 General safety instructions

/ DANGER

Please observe all warning information present on the device and in the instructions for use.

There is a danger to life if warnings are not obeyed. Severe personal injury or damage to property may occur.

Never use defective power cables

If a damaged or unsuitable power cable is used, it could result in a fire or electric shock. Only use power cables with PE contacts approved by the manufacturer.

Disconnect the power cable correctly

When disconnecting the power cable, always do so by holding the plug. Ensure that your hands are dry. There is a risk of electric shock.

Do not insert any objects into the housing

Objects inserted into the housing may result in an electric shock or damage to the device.

Do not place any objects on top of the device

If you place objects on top of the device, this can lead to overheating and fire.

Avoid penetration of liquid

Liquids seeping into the device may result in electric shock or device failure.

! CAUTION

Extensive damage to property may result if the device is not connected correctly

That is why you should observe the warning information:

Connection must be carried out by specialists

Please ensure that all steps are taken to avoid injuries or incorrect diagnoses.

- Only use the video cables specified by the manufacturer for the connection.
- Only use power cables with PE contacts.
- · Only use power outlets with PE contacts.
- Do not connect too many devices to a power outlet or extension cable.
- Observe the information provided by the respective manufacturer.
- If required by the application or local regulations, QA software must be used for quality control
 and documentation.

Connection in the USA and Canada

Molded power supply plugs must comply with the requirements for "hospital grade attachments" CSA Std. C22.2 No. 21 and UL 498.

Connection in China

Only use power cables approved for China. These power cables are identified by the labels "CCC" or "CQC".

Observe the country-specific regulations

Observe all regulations of the country in which the device is used.

NOTICE

Extensive damage to property may result if the device is not connected correctly

That is why you should observe the warning information:

- Desktop installation:
 - Place the device on a solid and level surface. The attached stand, as well as the installation surface, must be suitable for the weight of the device.
- For mounting on a wall or ceiling suspension:
 The mount unit must be suitable for the weight of the device.
- For installation in a rack:
 Observe the installation sequence, and provide ventilation for the device.

Provide adequate air circulation

When installing the device, ensure that there is adequate air circulation for operation. The permissible ambient temperature range must not be violated. Otherwise, the device could be destroyed by overheating.

Avoid sources of heat

Do not install the device in the vicinity of sources of heat, such as radiators, heating appliances or other devices that can generate or emit heat.

Do not subject the device to jolting or shocks

The device contains sensitive electronic components that could be damaged by jolting or shocks.

Only switch on a cold device following adaptation to room temperature

If the device is brought into a room with a higher or rising temperature, condensed water will form in and on the device. Do not switch on the device until the condensed water has evaporated. Otherwise, the device could be damaged.

2.1 General safety instructions

NOTICE

Extensive damage to property may result if the device is not connected correctly

That is why you should observe the warning information:

Transportation only in original packaging

Use the original packaging for transportation, and transport in the correct shipping position. Be sure in particular to protect the monitor LCD modules from shocks.

Care of device / cleaning agents

- Remove water drops immediately; extended contact with water discolors the surface.
- Only clean the surfaces using the cleaning agents referred to in the Instructions for Use.
- Monitor: The screen is extremely sensitive to mechanical damage. Absolutely avoid scratches, shocks, etc.

What to do if the device is faulty

If the following conditions exist, the device must be disconnected from the line power supply and checked by qualified personnel:

- Damage to the plug or power cable.
- After liquid seeps into the device.
- If the device has been exposed to moisture.
- If the device does not function or if a fault cannot be eliminated using the Instructions for Use.
- If the device has been dropped and/or the housing damaged.
- · If the device smells of burning or makes peculiar noises.

Be aware of the monitors aging

Note that monitors can fail as a result of aging, and that image properties such as brightness, contrast, and color value can change.

Do not touch the monitor screen

Due to mechanical pressure or electrostatic discharges, touching the screen can result in brief disturbances to the image.

2.2 Product-specific safety information

!CAUTION

Careful installation of the stand

The tilt, position, and height of the monitor can be changed when the monitor is installed on a stand. Note the following to prevent injuries or damage when installing or adjusting the stand:

- Make sure that you do not pinch or injure your fingers or any other parts of the body.
- Make sure the monitor does not collide with the table or other objects, which could result in damage.

/ CAUTION

Ensure monitor stability

Stability must be ensured after installing the monitor on a stand. An unstable stand can cause the monitor to tip over and cause injuries or damage. Therefore make sure the stand is stable.

- Only use a stand that has been tested for the monitor's weight and that is approved for tilting up to 10°.
- · The screw insertion depth into the monitor must be within the acceptable range.

NOTICE

Radio interference

This is a Class B device.

The device may cause radio interference or interfere with the operation of other devices in close proximity. In this case the user is encouraged to perform appropriate measures to correct the interference.

NOTICE

Medical System

- Do not connect devices that are not part of the medical system.
- If a stand is subsequently attached to the monitor, the stand must meet the requirements of standard EN 60601 so that the system, comprising the stand and monitor, meets the requirements.

Note

No zero error rate

LCD monitors do not have a zero error rate. For this reason, the image parameters can change over time, e.g. reduced luminance or changing/fading colors.

Note

Image quality

To maintain constant image quality, EIZO recommends cleaning the monitor on a regular basis and checking image properties in accordance with all applicable local regulations.

3 Description

3.1 Scope of delivery

The device and various components are included in the scope of delivery. After unpacking, check the scope of delivery for correctness and completeness.

Note

Keep the packaging material for subsequent transport of the device.

Device

The CuratOR LX491W is a Full HD 49" LCD Monitor. The monitor can be installed in a ceiling suspension, a wall mount, or a stand.

Product	Order number Description	
CuratOR LX491W (white)	6GF62602LA10	Full HD 49" LCD Monitor with white housing.
CuratOR LX491W (black)	6GF62602LA11	Full HD 49" LCD Monitor with black housing.

Components

The following components are included in the scope of delivery:

- 1 power cable for Europe
- 1 power cable for US
- · 1 power cable for China
- 1 power cable for Japan
- 1 DVI cable (3 m)
- 1 CD-ROM with the documentation
- 1 set of printed safety information
- 1 set of printed information on the battery per the German battery law BattG (in German)
- 1 remote control (AAA batteries included separately)

3.2 Monitor performance features

The CuratOR LX491W has the following features that permit a wide range of applications:

49" large screen diagonal

With a screen diagonal of 49" and a resolution of 1920 x 1080 pixels (2 MP), the LX491W is suitable for simultaneous use of several video sources. The monitor is particularly suited to the display of DICOM radiographic images, or as a secondary monitor in the surgical or endoscopy area.

LED backlight

The LX491W is equipped with an LED backlight optimized for bright environments. This enables a long service life even when operated at high luminance.

Fully Automated Stability

The LX491W has a Fully Automated Stability system that keeps luminance constant in accordance with medical standards such as DICOM or Gamma 2.2, for example. The integrated stability system ensures constant luminance using a built-in light sensor in the center of the backlight.

Preset Look Up Tables

The LX491W is precalibrated at the factory. A total of five practice oriented Look Up Tables (LUTs) have been preset. This calibration data makes installation and maintenance easier. As such, the monitor can be easily adapted to the respective application and local lighting conditions.

Adjustable LUT

The adjustable LUT functions allows manual recalculation of the lookup table. This allows users to modify the gamma model, luminance settings, or color parameters in accordance with local needs or preferences.

The LUT is recalculated via the OSD menu without having to access external calibration software. The reset LUT is calculated immediately, saved in the device, and used internally.

Force Mode

Using the Force Mode function, the LX491W can be adapted to special timing settings.

Video inputs

The LX491W can be connected to the imaging system using various video inputs, such as DVI, 3G-SDI, Composite, S-Video, and VGA. BNC connectors such as RGBS, SoG, and YPbPr can be inserted in the VGA input using an adapter cable. As a result, the LX491W can process digital and analog standard video signals (PAL, NTSC).

Any adjustments to the video settings that may be required can be carried out using an OSD menu (On Screen Display).

Simultaneous display of different video sources

The various video sources can be connected to the video inputs simultaneously and displayed on the monitor. The various signal inputs can be displayed either as a "Picture in Picture" (PiP) representation, or next to one another as "Picture and Picture" (PaP) using the widescreen format.

Protective screen

The LX491W has an anti-glare protective screen fitted in front of the LCD panel to protect the surface from bumps and scratches. The monitor with protective screen is protected at the front against moisture (IP65 degree of protection).

The space between the protective screen and the panel is sealed to prevent dust from entering, thus helping ensure the internal surfaces remain clean.

4 Setup and installation

/!\CAUTION

Changes to device

Do not make any mechanical or electric changes to the device. Otherwise the device warranty becomes invalid.

The manufacturer is not liable for changes made to the device.

4.1 Installation site

NOTICE

The power switch and connections must be accessible at all times

When installing and connecting the monitor, ensure that the power switch and the connections are accessible at all times.

NOTICE

Condensation

If the device is brought into a warm environment from a cold one, condensation may form in the device. This could result in a short circuit when switching on the device, damaging it.

 Wait until the condensed water has evaporated, including that inside the device, before you switch it on. This can take several hours.

NOTICE

Overheating

Ventilation holes are located on the rear of the housing.

If the ventilation holes are covered or closed, the heat generated in the monitor will not be dissipated sufficiently.

- · Do not cover the ventilation holes.
- Do not close the ventilation holes.
- The minimum distance from the back and side of the monitor to the wall must be 10 cm, and at least 15 cm from other devices.
- The ambient temperature must be in the acceptable range of +5 °C to +40 °C.

NOTICE

Dusty environment

The monitor is intended for use in the clean environment of medical diagnostics. In dusty environments, ventilation holes in the back can allow dust to penetrate into the monitor.

In the worst case, deposits are possible which become evident as dark spots in a white picture and result in deterioration of the luminance.

 Protect the monitor from dust, for example through appropriate construction measures at the installation site.

Note

Reflections on the screen

The monitor has an anti-glare surface that is only effective if the screen is clean and grease-free.

- · Comply with the specifications for cleaning.
- Position the monitor to avoid reflections on the display area.
 Reflections can be caused by lights, windows, furniture with shiny surfaces, or light-colored walls.
- In order to reduce reflections on the monitor, only use non-dazzling reflector bulbs for the ceiling lighting.

Note

Shocks and impacts

The monitor is sensitive to mechanical influences. Shocks and impacts on the panel surface can lead to total failure.

· Ensure that such mechanical influences at the installation site are avoided.

Note

Movable installation

If the monitor is installed such that it can move, make sure that persons or objects in the facility are not endangered by the monitor's range of movement.

Note

During transport, use the original packaging or service packaging.

4.2 Installing the monitor

NOTICE

Holders

- Mounts must be tested and approved by the manufacturer for the weight to be supported.
- An installed stand must be sufficiently stable to withstand tilting of up to 10° without toppling the monitor.

Note

Inaccessible cover

Before installing the monitor in the selected holder, check whether the screws for removing the cover of the connection panel are accessible. If not, remove the connection panel cover prior to installation.

The monitor has a VESA 400x200 mounting interface and can be installed in a suitable ceiling suspension or wall mount.

Note the following during installation:

- The maximum torque for attaching to the holder is 10 Nm.
- The screws used for attaching to the holder must meet the following requirements:

No.	4		
Thread	M8		
Strength	8.8 in accordance with ISO 898-1		
Insertion depth	10 mm to 30 mm		

5 Connecting

5.1 Safety information for connection

All safety information and warnings for the device must be observed to ensure danger-free operation.

/ CAUTION

Changes to device

Do not make any mechanical or electric changes to the device.

EIZO GmbH is not liable for changes made to the device.

!CAUTION

Shielding measures

Follow all shielding measures in accordance with local EMC directives. If these guidelines are not observed, device malfunction may result.

!CAUTION

Grounding

The permissible leakage current is not exceeded during the first fault event in accordance with EN60601-1. The device is grounded with an additional protective conductor to ensure the greatest possible electric safety.

/ CAUTION

Excessive currents, short circuits, and ground faults

In accordance with national standards and regulations, protection against excessive currents, short circuits, and ground faults must be incorporated into the building installation.

NOTICE

Changes to device settings

Device settings may only be adjusted by service personnel.

NOTICE

Disconnecting from line power

Always set the power switch to "Off" before disconnecting the device from power. Otherwise the device could be destroyed.

NOTICE

Cable installation

Observe the following instructions:

- · Only shielded cables are to be used for all signal connections.
- If the relevant facility is available on the connector, all plug connections must be screwed tight or locked.
- · The connecting cables must not be kinked.
- The minimum bending radius of a connecting cable generally equals five times the cable diameter.
- Do not route signal cables and power cables next to one another. Otherwise, line power subject to heavy interference could result in reversible pixel errors.
- The device must not share a line power supply with motors or valves (interference!).
- Externally connected cables can represent a trip hazard. Make sure that all incoming cables are safely routed.
- If the device offers strain relief mechanisms for the cables, use them to prevent unintended loosening of connected cables.

5.2 Connector locations

The connectors are located in the connection panel behind a cover on the back of the monitor. The power switch is not covered and is freely accessible.

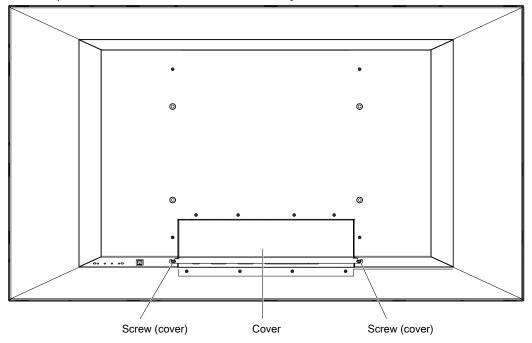


Fig.: Rear view with cover

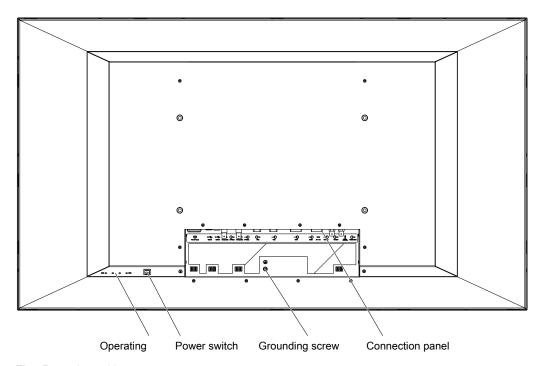


Fig.: Rear view without cover

Note

An additional protective conductor may be connected to the grounding screw.

5.3 Connection panel

The connection panel is visible after removing the cover. The LX491W has connectors for power, communication, and video input and output signals.

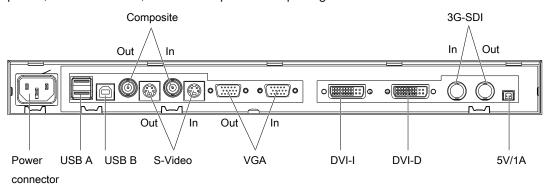


Fig.: LX491W connection panel

5.3 Connection panel

Signal inputs

The following signal inputs are available on the connection panel.

- VGA input (15-pin D-sub socket): In addition to a standard VGA signal, other analog signals can be connected with a suitable adapter. Examples include component, composite, and RGB video using a multi-BNC to VGA adapter.
- DVI-I and DVI-D input
- 3G-SDI (BNC)
- S video with Y/C signal (4-pin mini-DIN socket)
- Composite (BNC)

Note

Interlaced signals

If an interlaced signal is applied to a signal input (except for the 3G-SDI input) and the "Deinterlacing" function in the "Signal" OSD menu is set to "Normal", brightness is reduced by approx. 15%.

The "Field to frame" and "Static mesh" deinterlacing settings cause a slight reduction in brightness only, but can lead to comb effects.

Also refer to the "Deinterlacing" function in the "Signal" menu [▶ 42].

Signal outputs

The following signal outputs are available on the connection panel.

- VGA output (15-pin D-sub socket):
- 3G-SDI (BNC)
- S video with Y/C signal (4-pin mini-DIN socket)
- Composite (BNC)

Additional connectors

- Power connector (appliance plug)
- USB: 2 downstream (type A) and 1 upstream (type B) connector.
- Voltage supply: 5V/1A connector to supply external devices.

5.4 Connecting the signal cable

!CAUTION

Opening the connector panel cover

Only service may open the connector panel cover. Patients must not be present when the cover is open.

ACAUTION

Connector

Connectors may only be plugged in or removed by Service when the device is switched off.

NOTICE

Cable

- The picture quality, interference immunity, and emitted interference of the entire system depend on the cable quality and length.
- When using a BNC to VGA adapter cable, the signal cables, for example, red, green, and blue, must be of the same length to prevent a loss in image quality.
- · Use only the cables specified by EIZO or the transmission links available from EIZO.

Note

Video source settings

- The video source is set on the monitor using the EDID data transmitted via the DDC interface. If the video source cannot interpret the EDID data, the monitor automatically attempts to adjust to the signal clocking of the video source.
- Do not change these settings. Otherwise, the images will not be displayed correctly.

The connection panel for the signal cables and power supply is located on the back of the monitor. All signal inputs can be connected simultaneously.

Prerequisite

- The monitor is installed in the ceiling suspension or wall mount unit, or on a stand.
- The screws for attaching the cover to the connection panel are accessible.
 See also Connector locations [▶ 18] and Installing the monitor [▶ 16].

Procedure

- 1. Using a suitable tool, remove the connection panel cover.
- 2. Connect the signal cables to the monitor.

Note: All signal inputs can be connected simultaneously.

- 3. Tighten the screws to secure the signal cables.
- 4. Place the cover back on the connection panel.

5.5 Connecting the power cable

5.5 Connecting the power cable

! DANGER

Connecting to line power

- The device is designed for line power with a grounded neutral conductor.
- To avoid risk of electric shock, this device must only be connected to line power with a
 protective conductor.
- Contact the responsible building technician or a qualified electrician if you are uncertain whether the line power is equipped with a protective conductor.

!CAUTION

- Only use power cables or device connection cables with protective conductors and appliance plugs according to DIN 49 547, IEC 60320 (max. length 3 m). Furthermore, the cable must adhere to all local safety regulations applicable to the specific country in which the monitor is used.
- Device fuses should only be replaced in repair centers or by service personnel.
- Note for North America: Molded power plugs must comply with the requirements for hospitals with respect to CSA Std. C22.2 No. 21 and UL 498.

/ CAUTION

Opening the connector panel cover

Only service may open the connector panel cover. Patients must not be present when the cover is open.

/ CAUTION

Connector

Connectors may only be plugged in or removed by Service when the device is switched off

The monitor is supplied power via an appliance plug in the connection panel on the back of the monitor.

Prerequisite

- The monitor is installed in the ceiling suspension or wall mount unit, or on a stand.
- The screws for attaching the cover to the connection panel are accessible.
 See also Connector locations [▶ 18] and Installing the monitor [▶ 16].

Procedure

- 1. Using a suitable tool, remove the connection panel cover.
- 2. Connect the appliance plug to the monitor power socket.
- 3. The power cable can be secured using a cable grip.
- 4. Place the cover back on the connection panel.

6 Commissioning

Note

Factory settings

All monitors are optimally preset in the factory, such that changes are not usually required.

6.1 Switching on the monitor and video source

Note

To obtain the best possible results, the video source should support communication via the Display Data Channel (DDC).

The monitor and video source can be switched on in any order.

Switch on

- · Switch on the monitor.
 - The EIZO logo appears briefly.
- Switch on the video source.
 - If there is no signal the monitor displays a black screen.
 - If there is a signal, the image is displayed accordingly.

Note

Image not displayed?

If an image is not displayed after the monitor is switched on and a video signal is present:

 check the system for basic connection and operating errors before contacting service personnel.

6.2 Avoiding image sticking

Image sticking may occur with LCD monitors. Image sticking is an effect whereby a faint image of the previous screen contents can be seen after the display contents have changed.

The following measures can reduce or prevent image sticking:

- Use a screen saver with regularly changing images
- Switch off the monitor when it is no longer needed.
- The monitor has an energy saving mode:

 If the application in use supports the energy saving mode, activate it.

6.3 Check for pixel defects

Note

Energy saving (Power Management)

The monitor supports various energy saving settings, called Power Management (PM). When PM is active, the monitor backlight switches off automatically for example, if the monitor is without a video signal for an extended period.

Also observe the operating system manufacturer's instructions regarding power management settings.

6.3 Check for pixel defects

Pixel defects (small bright or dark dots) can occur in LCD monitors. During the manufacturing process, all monitors are checked for the permitted number of defective pixels.

Defective pixels cannot be corrected.

6.4 Setting the image geometry

The monitor automatically recognizes the standard used, and has preprogrammed set-up values for each standard. However, depending on the graphics card used, it may still be necessary to align and size the picture for the selected standard.

6.5 Adjustment of monitor and video source

Adjustment of resolution and refresh rate

Every LCD monitor has certain limits, such as maximum resolution and refresh rate.

 Set the graphics card for operation of the monitor in such a way as to maintain the limit values.

Fine adjustment of brightness and contrast

/!\CAUTION

The exact brightness and contrast settings are only possible with a photometer.

/!\CAUTION

Fine adjustment of analog inputs: Only via 15-pin D-sub and DVI-I sockets Fine adjustment for digital input: not necessary

- Perform the fine adjustment of the monitor only via the two analog ports (15-pin D-sub and analog signals on the DVI-I socket).
- Digital signals on the DVI-I input do not require a fine adjustment of brightness and contrast since the signal is always displayed in an optimal manner. With a fine adjustment, it is possible that gray scales are not displayed.

RGB video sources (via 15-pin D-sub or DVI-I connector) supply analog signals that are basically intended for conventional CRT monitors and that are processed directly by them.

In contrast, the analog signals must be converted into digital signals for the LCD monitor by a video digitizer. Depending on the video source, cable length, and video mode (e.g. VGA, SVGA, XGA), this conversion may cause certain deviations that cannot be fully corrected by the monitor automatically.

 Perform a manual fine adjustment to obtain the optimum image display with an analog input signal (VGA/DVI-I). During the manual fine adjustment, fine tune the monitor (or video digitizer to be more precise) to the respective video source.

Adjusting the monitor for optimum performance with a graphics card

In order to optimize the monitor settings for the installed graphics card and guarantee that all gray scales can be distinguished, we recommend only adjusting the brightness and contrast for the analog inputs.

Please note that these settings have no impact on the calibration in the look up table: EIZO GmbH monitors are calibrated at the factory and maintain these settings:

• To reduce the brightness using the OSD operator controls, use a picture with 0 % gray scale (black) and a suitable measuring instrument.

Reduce the brightness until the measuring instrument displays constant values: the measured value no longer changes.

Then increase the brightness slightly until the monitor is just above the lowest black level (one step is generally sufficient).

 To set the white value, use a test pattern with 100 % gray value (white) and a suitable measuring instrument.

To ensure the black value remains unchanged, set the contrast only.

 Increase the contrast until the measuring instrument no longer registers an increase in luminance.

Then adjust the contrast regulator setting 1 or 2 steps back until the measured brightness is just under the maximum value.

• Make sure once again that the black value has not changed.

The black value can change if significant corrections were made to the contrast. In this case, repeat the steps above.

Result

The monitor is now configured for optimum performance with the installed graphics card. If the luminance is still not satisfactory, you can increase the black and white values further by adjusting the backlight in the OSD menu.



Permanently higher setting for backlight

The permanently higher setting for backlight reduces the service life and brightness performance.

6.6 Adapting the monitor using Force Mode

Setting without a measuring instrument

Precise adjustment is only possible with a measuring instrument.

If there is no measuring instrument available for fine adjustment, proceed as follows:

- · Use the SMPTE test pattern.
- Adjust the brightness so that image sections with 5 % and 0 % gray value still visibly contrast.
- Adjust the contrast so that image sections with 95 % and 100 % gray value still visibly contrast. To adapt the luminosity to the ambient lighting, adjust the backlight luminance.
 CAUTION: the factory-set brightness is no longer maintained.

Note

To set the image parameters you can also use the adjustable LUT functions in the "Adjust LUT" menu [> 51].

Automatically setting the video source

In order to automatically adjust the video source, it must support communication via the Display Data Channel (DDC), and the devices must be connected correctly. When switched on, the monitor's EDID data (Extended Display Identification Data) is read out, and the video source can detect the monitor.



Installation and parameterization of the video source

Please refer to the video source manufacturer's manual for detailed information about installation and configuration of the video source.

6.6 Adapting the monitor using Force Mode

NOTICE

Force Mode is a service tool

Force Mode should only be used by service personnel to determine an unknown timing and implement it in the monitor.

6.6.1 Introduction to Force Mode

Various factory-set timings are saved in the monitor. As soon as a video signal is connected, an appropriate timing is sought. "Auto In Process" is displayed during this phase. These timings are compatible with the standard video signals provided by current graphics cards.

If no image or only an unclear image is displayed on the monitor, the signal is outside the standard supported range. Such signals are common with older medical equipment. A timing of this sort can be set using the Force Mode function.

In many cases, a signal outside the standard range is displayed correctly without having to use Force Mode. This is possible because a large number of known timings are saved in the monitor.

6.6.2 Basic information on timing

Scanning in the interlaced procedure and procedures without line interlacing

There are two different scanning systems. They differ in the technology used to display the image on the screen TV signals and monitors that are compatible with them are normally set to the interlaced procedure; computer signals and monitors compatible with them are normally set to the procedures without line interlacing (progressive or non-interlaced) These two formats are not compatible; one of them must first be converted before the signals can be processed together.

In the case of interlaced scanning, each image is divided into two separate fields. An image therefore comprises two fields. An interlaced image is output on the screen in two scans. The horizontal lines of the first field are scanned first, and then, again starting at the top of the image, the horizontal lines of the second field are scanned between the first set of lines. Field 1 consists of the lines 1 to 262 1/2, and field 2 of the lines 262 1/2 to 525. Only a few lines are displayed at the top and bottom of each field.

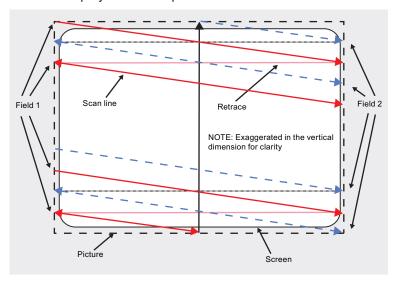


Fig.: Interlaced scanning system

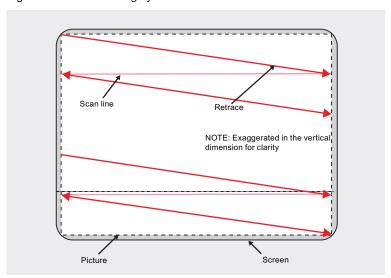


Fig.: Progressive (non-interlaced) scanning system

6.6 Adapting the monitor using Force Mode

A progressive image is output on the screen in that all horizontal lines are scanned from top to bottom in one scan.

Horizontal timing diagram (the vertical timing diagram is identical)

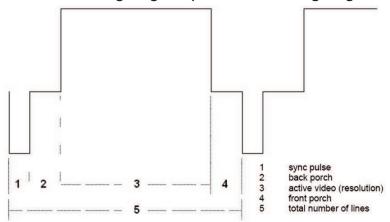


Fig.: Timing diagram


Force Mode menu	Timing diagram		
Horizontal Resolution	3 - Active video (horizontal resolution)		
Vertical Resolution	3 - Active video (vertical resolution)		
Total horizontal lines	5 - Total number of lines (horizontal)		
Horizontal blank pixels	4 – Front Porch (horizontal)		
Vertical blank pixels	4 – Front Porch (vertical)		

6.6.3 Entering timing data

If you know or have already determined the timing data, you can enter it directly in the OSD "Force Mode" menu.

Entering known timing data

- 1. Open a test pattern with a clearly defined border (e.g. SMPTE image).
- 2. Open the "Force Mode" menu is the OSD.

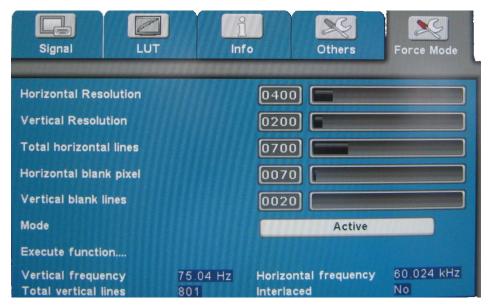

- 3. Use the keys to enter the timing data in the input fields.
- 4. Set the "Mode" to "Active".
- 5. Apply the entered data by selecting "Execute function..."

Note

Always save the data with "Execute function...".

The timing data are only entered into memory with the "Execute function..." command, and the image is then displayed with the new parameters.

6.6.4 Determining the timing data

If you know none or only part of the timing data, you can determine the data incrementally in the "Force Mode" menu.

Note

All settings must be implemented with the zoom factor "1 to 1" ("Geometry" menu).

Sequence for determining a timing

- 1. Determine the start values using the auto function.
- 2. Enter the known data.
- 3. Optimize the scanning frequency.
- 4. Optimize the geometry and resolution.

Step-by-step instructions for setting a timing can be found in the next sections.

Determine start values using auto function

The following commands can be used to trigger automatic determination of the timing data, and provides rough approximation of the applied timing.

- 1. Open a test pattern with a clearly defined border (e.g., SMPTE image).
- 2. Open the "Force Mode" menu in the OSD.



- 3. Set "Mode" to "Start".
- 4. The start values are determined when you select "Execute function..."

 Note: The timing data are only entered into memory with the "Execute function" command, and the image is then displayed with the new parameters.

After the "Execute function..." "Mode" switches to "Active" and you can optimize the start values.

Entering known data


If some of the timing data are known they can be entered after determining the start values in order to simplify the subsequent steps.


If the resolution is known for example, it can be used for the approximation of "Total horizontal lines". In this case, "Total horizontal lines" has to be greater than "Horizontal Resolution". If "Total horizontal lines" is less than "Horizontal Resolution" during setting, that latter has to be reduced. A larger "Horizontal Resolution" can result in the image being split vertically.

Optimize scanning frequency

The correct signal frequency now has to be set with "Total horizontal lines".

- 1. To set "Total horizontal lines" correctly, open a test pattern with a clearly defined border (e.g., SMPTE image).
- 2. Correct the signal frequency so that no interference is present in the picture. If the distance between interferences becomes larger, and these therefore also become fewer, you are changing the values in the correct direction.

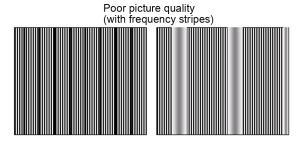


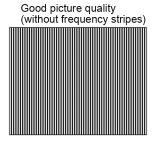

Optimize geometry/resolution

The active resolution range is defined by the values of the "Horizontal Resolution" and "Vertical Resolution" input fields. This range is displayed in the center of the monitor in black.

The inactive resolution range is the unused range between the active resolution and the maximum resolution of the monitor. In Force Mode, this range is displayed gray (gray scale monitors) or blue (color monitors).

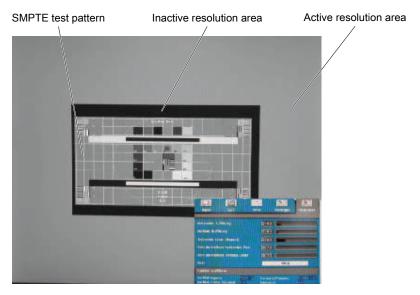

To adapt the active resolution range to the actual video signal resolution, perform the following steps:

- 1. Open a test pattern with a clearly defined border (e.g., SMPTE image).
- 2. In the "Force Mode" menu, use the "Horizontal blank pixel" input field to shift the left edge of the SMPTE image pixel-exact to the left inside edge of the gray/blue area.
- 3. Correct the values in the "Horizontal Resolution" input field until the right internal border of the gray/blue area is pixel-exact at the right edge of the SMPTE image.

6.6 Adapting the monitor using Force Mode

- 4. Use the "Vertical blank lines" input field to shift the top edge of the SMPTE image pixel-exact to the top internal border of the gray/blue area.
- 5. Correct the values in the "Vertical Resolution" input field until the bottom internal border of the gray/blue area is pixel-exact at the bottom edge of the SMPTE image.
- 6. Apply the entered data by selecting "Execute function..."

Note: Gray/blue range not visible

If no gray/blue area is visible on the left and right sides of the image, reduce the "Horizontal Resolution" until the gray/blue area becomes visible. If no gray/blue area is visible on the top and bottom sides of the image, reduce the "Vertical Resolution" until the gray/blue area becomes visible.

Timing successfully set

The timing was set successfully. You can now exit the OSD menu. The monitor will recognize the newly set timing and select it whenever connected.

Note

Only one timing can be saved in Force Mode itself.

Note

Once all settings in Force Mode have been performed, the zoom factor can be changed in the "Geometry" menu.

7 Operation

This section describes the operator controls and functions of the monitor OSD menu.

!CAUTION

Settings in the OSD menu

Changes to the OSD menu settings can impact the image properties. For this reason, OSD menu settings should only be changed by trained personnel.

Measures in the event of a failure


Note

Device malfunction in operation

If the device is not working properly, check the system for basic connection and operating errors before contacting service personnel.

7.1 Operator controls

The CuratOR LX491W has four keys for OSD operation located in the lower right corner.

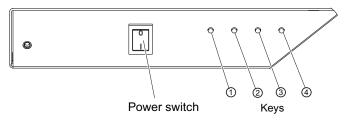


Key functions

The keys have the following functions within the OSD menu:

Key	Function					
1	Open OSD menu					
	Toggle					
2	Navigate between the submenu / tabs					
	Decrease/change value					
3	Navigate between the submenu / tabs					
	Increase/change value					
4	Exit OSD or submenu (settings are retained)					
Power switch	Switch the device on and off					

7.2 Remote control

The CuratOR LX491W is delivered with a remote control.

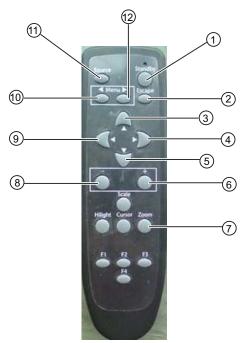

Standby
 Down
 Escape
 Plus
 Menu L
 Up
 Zoom
 Source
 Right
 Minus
 Menu R

Key	Action				
Menu L	Open the OSD menu Novigate tab to the left				
Marris	Navigate tab to the left				
Menu R	Open the OSD menuNavigate tab to the right				
Escape	Exit OSD				
Up	Navigate upwards in menu				
Left	Close submenu				
Right	Open submenu				
Down	Navigate downwards in menu				
Minus Decrease value					
Plus	Increase value				
Standby	Switch the device to standby mode (this key must be held down for an extended period to avoid inadvertently turning off the monitor when accidentally pressing the standby button)				
Source	Switch between connected signal sources				
Zoom Change the zoom ("One to One", "Fill all", "Fill to aspect ratio")					

7.3 Lock or unlock OSD menu


!CAUTION

Locking and unlocking the OSD menu

- Only authorized service personnel may lock or unlock the OSD menu.
- The OSD must be locked if inappropriate operation by the user can impact the intended use of the monitor.

The OSD menu of the LX491W is locked when shipped from the factory.

To lock or unlock the OSD menu, proceed as follows:

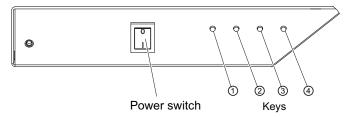

- 1. Press key @ once.
- 2. Then press key 2 three times.

The OSD menu is now locked or unlocked, depending on its initial state.

7.4 Picture layout (PaP, PiP, PoP)

PaP: "Picture and Picture" arrangement

Two images from a primary and secondary video source are displayed next to one another.

PiP "Picture in Picture" arrangement

The image content of the secondary video source is displayed on the main image, the primary video source.

PoP "Picture over Picture" arrangement

NOTICE

PoP image arrangement not for live operation

When working with the PoP image arrangement, there is a time delay when displaying the images (snapshots).

• Do not use the PoP image arrangement for live operation.

Up to five images of various video sources may be displayed next to one another.

- The image from the primary video source is displayed in the larger segment in the lower section of the window.
- The other images are displayed in a smaller format in the upper section of the window.

7.4 Picture layout (PaP, PiP, PoP)

The primary video source can be selected. The four additional video sources are added automatically in accordance with the following table of compatible video source combinations:

PiP / PaP	DVI 1	DVI 1	DVI 2	VGA	S-Video	SDI	Composite Video
Compatible combinations	digital	analog	digital				(CVBS) / SoG
DVI 1 digital	-	-	0	-	0	01	0
DVI 1 analog	-	-	0	-	0	01	0
DVI 2 digital	0	0	-	0	-	-	-
VGA	-	-	0	-	0	0	0
S-Video	0	0	-	0	-	-	-
SDI	01	01	-	0	-	-	-
Composite Video (CVBS) / SoG	0	0	-	0	-	-	-

Video sources can be displayed simultaneously

Image settings on the PiP/PaP or PoP window

The image display and quality of the various PiP / PaP / PoP windows can be set individually.

- When "Picture layout PaP; PiP; PoP" is enabled in the "Signal" menu, an "Active adjustment window" selection appears in the main OSD menus for "Picture", "Image", and "Signal".
- The settings for the OSD window can be applied to the main or secondary window.

Note

Picture layout PaP, PiP, PoP

- In order to adjust the PiP or PaP window settings in the "Image" menu, "Active adjustment window" must be set to "PiP/PaP Window (2)".
- The number displayed on the right side of the OSD window indicates which window is currently active for the setting ("1": main window, "2": secondary window)

⁻ Video sources cannot be displayed simultaneously

¹ **Note:** Flickering horizontal stripes may occur when the OSD is open.

7.5 Overview of OSD menu

You can use the OSD menu to change monitor settings or display information about the monitor.

If no input signal has been applied, the OSD only offers a limited set of functions.


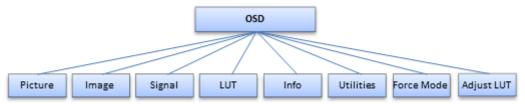


Fig.: Overview of the OSD menu

7.6 "Picture" menu

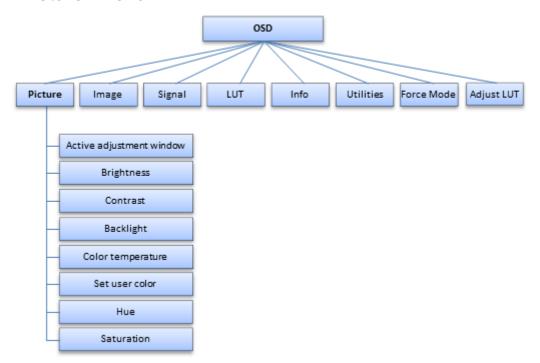


Function	Values	Description
Active adjustment	Main window (1)	Selecting the active adjustment window
window	PiP window (2) Default: Main window (1)	Note: The function is displayed only if the "Picture layout PaP, PiP, PoP" function is active within the "Signal" menu.
		The "Active adjustment window" establishes the window for which the settings in this menu apply.

Operation 7.6 "Picture" menu

Function	Values	Description
Brightness	Slider control	Set brightness
	Default: 50	Adapting the representation of darker picture areas.
		Note: The brightness settings are already optimized for digital DVI signals. Do not change these values manually, as this can negatively impact picture quality (loss of gray scales).
Contrast	Slider control	Set contrast
	Default: 50	Adapting the representation of brighter picture areas.
		Note: The contrast settings are already optimized for digital DVI signals. Do not change these values manually, as this can negatively impact picture quality (loss of gray scales).
Backlight	Slider control	Adjust brightness of monitor backlight
		Adjustment of overall brightness to ambient lighting.
		Note: The backlight settings always apply for all windows.
Color temperature	Native	Set the desired color temperature or hue
	9300°K 7300°K	The color temperatures can be selected from among three fixed settings or one adjustable setting.
	User Default: Native	It is possible to define a different color temperature for each video input.
Set user color		Define variable color temperature
		The color setting defined here can be set in the "Color temperature" menu item using the "User" setting. Note: If the color coordinate is adjusted with this function, some color levels may be lost. To prevent a loss in color levels, calibrate the LUT color point with suitable QA software.
Red	Slider control Default: 50	Select red component of display.
Green	Slider control Default: 50	Select green component of display.
Blue	Slider control Default: 50	Select blue component of display.
Hue	Slider control Default: 50	Set the color hues for RGB and video signals. Note: "Hue" may only be set for S-Video or Composite signals.
Saturation	Slider control Default: 50	Set the color saturation for RGB and video signals. Note: This menu is not displayed for the monochrome CVBS signal.

7.7 "Image" menu

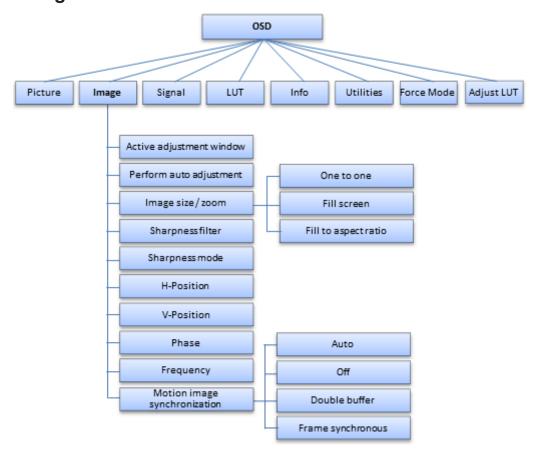


Function	Values	Description
Active adjustment win-	Main window (1)	Selecting the active adjustment window
dow	PiP/PaP window (2)	Note: The function is displayed only if the
	Default: Main window (1)	"Picture layout PaP, PiP, PoP" function is active within the "Signal" menu.
		The "Active adjustment window" establishes the window for which the settings in this menu apply.
Perform auto adjust-	Execute	Perform auto adjustment
ment		"Position", "Phase", and "Frequency" are set automatically upon confirmation of the selected function. Note: "Perform auto adjustment" can only be used for analog signal inputs.

Operation 7.7 "Image" menu

Function	Values	Description
Image size / zoom	1:1 Fill to aspect ratio Fill screen Default setting: 1:1	Selection between different picture size settings: One to One: The picture is displayed on screen in its original resolution. Fill to aspect ratio: The picture is zoomed to the maximum screen area with retention of the aspect ratio. Fill screen: The picture is displayed to fill the complete screen (1920 x 1080 pixels).
Sharpness filter	Slider control	Select the sharpness setting Note: To display the "Sharpness filter" function, "Sharpness mode" must be set to "Normal". Use the sharpness setting to reduce scaling artifacts or make the image "softer". A lower value creates a "soft" image impression and a higher value a "sharp" one. Note: The optimum sharpness setting has to be determined visually. Common filters are available for the RGB sources (VGA, DVI). The interpolation filters depend on the input resolution. At lower resolutions, the filter calculates the value for the non-controlled pixels.
Sharpness mode	Enhanced Normal Default setting: Enhanced	Setting image sharpness mode Select "Enhanced" to display the image with the sharpest setting. Under "Normal", image sharpness can be set using the "Image sharpness filter" slider control.
H-position	Slider control	Shift picture in horizontal direction
V-position	Slider control	Shift picture in vertical direction
Phase Frequency	Slider control Slider control	Set the frequency and phase of the input signal Correct potential blur in the vertical lines by adjusting the "Frequency/Phase" setting. Note: Use a vertical line in a "Pixel On/Off" test pattern to make the adjustment.

Function	Values	Description
Motion image synchro-	Auto	Optimize video display
nization	Off Double buffer	Using this function, timing-based video artifacts (tearing) can be avoided.
	Frame synchronous Default setting: Auto	Note: This function is displayed only if the "Picture in Picture (PiP)" setting is selected for "Layout format" in the "Signal" menu. Only settings that can be shown technically given the current use of the display can be selected. The function setting always applies to the active setting window
		Auto: The best possible automatic adaptation to the applied signals.
		Off: Motion image synchronization is switched off and the monitor is running in asynchronous mode. Tearing is not pre- vented. Lowest latency between input and output signal.
		Double buffer: Tearing is prevented. The latency between input and output is 1-2 frames. Recommended setting for signal frequencies greater than 60 Hz or less than 50 Hz".
		Frame synchronous: Tearing is prevented. Latency less than 1 frame. Recommended setting for signal frequencies of 50 to 60 Hz ⁻¹ .
		") The "Picture layout PaP, PiP, PoP" function must be disabled.

7.8 "Signal" menu

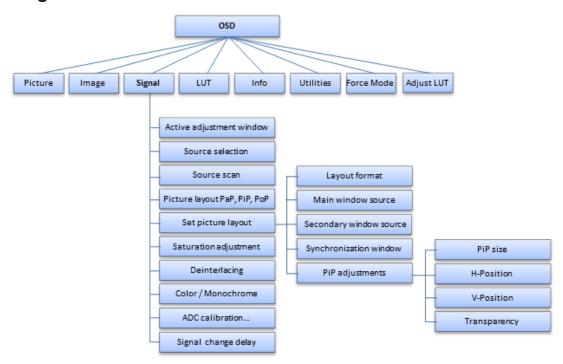

Function	Values	Description
Active adjustment	Main window (1)	Selecting the active adjustment window
window	PiP window (2) Default: Main window (1)	Note: The function is displayed only if the "Picture layout PaP, PiP, PoP" function is active.
		The "Active adjustment window" establishes the window for which the settings in this menu apply.
Source selection	VGA / DVI 1 Analog /	Select source for main display
	DVI 1 Digital / DVI 2 Digital / CVBS / S-Video / SDI /	Note: This function is only displayed when the "Source scan" function is set to "Disabled".
		Selection of picture source for full format image.
Source scan	Enabled / Disabled	Activate / deactivate automatic source scan
	Default: Enabled	Enabled: If the displayed source is no longer available, the monitor searches au- tomatically for the next available source.
		Disabled: If the displayed source is no longer available, no image is displayed: The monitor does not search for available sources. The monitor screen appears black. An image is not displayed again unless a signal is applied to the set source or a different signal source is selected manually.


Function	Values	Description
Picture layout PaP, PiP, PoP	Enabled / Disabled Default: Disabled	 Switch the preconfigured picture layout (PaP; PiP; PoP) on or off When "Picture layout PaP, PiP, PoP" is enabled, an "Active adjustment window" selection appears in the main OSD menus for "Picture", "Image", and "Signal". The settings for the OSD window can be applied to the main or secondary window. Note: The number displayed on the right side of the OSD window indicates which window is currently active for the setting ("1": main window, "2": secondary window)
Set picture layout		
Layout format	Picture in Picture (PiP) Picture and Picture (PaP) Picture over Picture (PoP) Default: Picture in Picture (PiP)	The following picture layouts can be selected using this function: Picture in Picture (PiP): The image content of the secondary video source is displayed on the main image (primary video source). Picture and Picture (PaP): Two images from a primary and a secondary video source are displayed next to one another. Picture over Picture (PoP): Up to five images of various video sources may be displayed next to one another. The image from the primary video source is magnified and displayed in the lower section of the window. The other images are displayed in a smaller format in the upper section of the window. Note: For "Picture over Picture (PoP) images" there is a time delay when displaying images. This setting may not be used for live mode.
Main window source		Select the video sources displayed in the respective "Picture in Picture (PiP)", "Picture and Picture (PaP)", and "Picture over Picture
2nd Window source		(PoP)" windows. See also Picture layout (PaP, PiP, PoP) [▶ 35] for compatible video source combinations.
Synchronization window	Main window (1) PiP/PaP window (2) Default: Main window (1)	Select whether the main or secondary window serves as the synchronization window. Note: The "Synchronization Window" menu is only available when the "Picture in Picture (PiP)" setting is active.

Operation 7.8 "Signal" menu

Function	Values	Description
PiP adjustments	The "PiP adjustments" menu only appears when the "Picture in Picture (PiP)" setting is active. Note: When initially selecting sources the "PiP adjustments" function is not always displayed. In this case, change to another OSD menu (such as "LUT") and then back to the "Signal" menu. The "PiP adjustments" function is now displayed.	
PiP size		Select the "Picture in Picture (PiP)" window size
H-position	Slider control	Set the horizontal position of the "Picture in Picture (PiP)" window
V-position	Slider control	Set the vertical position of the "Picture in Picture (PiP)" window
Transparency	Slider control	Select the "Picture in Picture (PiP)" window background ("opaque" or "transparent")
Saturation adjustment	Enabled / Disabled Default: Disabled	Activates the color adjustment of "Saturation" in the main "Picture" menu.
Deinterlacing	Field to frame	Setting the "Deinterlacing" method
	Static mesh Normal Default: Normal	Note: The "Deinterlacing" function is shown only in conjunction with a supporting signal source. Note: If an interlaced signal is connected to a signal input (except for the 3G-SDI input) and the "Normal" setting is selected, brightness is reduced by approx. 15% due to the color space conversion of the deinterlacing. The deinterlacing settings "Field to frame" and
		"Static mesh" use other deinterlacing methods and cause only a slight reduction in brightness. The "Static mesh" setting is optimized for still images, but could cause comb effects in moving images. The "Field to frame" setting is optimized for moving images, but could cause comb effects in still images.
Color / Monochrome	Color Monochrome YPbPr/YCbCr Default: Color	Switch signal between monochrome and color If "Color" is set, a monochrome signal is displayed green. For correct gray scale representation, select the "Monochrome" setting.


Function	Values	Description
ADC calibration	Execute	Automatically calibrate A/D converter for the applied video level
		Note: This function is only displayed when the "Saturation adjustment" function is disabled.
		The video level range of the system is checked, and the monitor set accordingly. This results in optimum adjustment of the individual RGB A/D converters to the video source.
		The calibration results in a larger and more flexible video level range (e.g., the full brightness is also achieved at 700 mV if the video level is limited at the value of 700 mV).
		Note: A suitable test pattern is required. The A/D converters have been factory-set and need not be readjusted. The function is shown only in conjunction with a supporting signal source.
Signal change delay	Slider control	Note:
	Default: 1	The "Signal change delay" menu item is not available unless "Source scan" is set to "Enabled".
		"Signal change delay" can only be set when "Picture layout PaP, PiP, PoP" is set to "Disabled".
		The following four parameters are relevant for the function:
		H frequency
		V frequency
		V total
		Interlaced/non Interlaced
		If one of these parameters changes, the monitor treats it as a timing change and initiates resynchronization via an "Auto in progress". To prevent this from happening as a result of each and every minor signal disturbance, the value representing the permissible number of faulty or changed borders must be increased in the case of unstable signal sources.
		Note: The higher the slider setting, the longer it will take for a desired timing change to occur (delayed by a number of milliseconds).

7.9 "LUT" menu

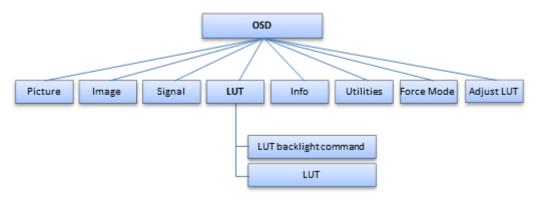


Function	Values	Description
LUT backlight command	Enabled	Control backlight using LUT
	Disabled Default: Enabled	When "Enabled" is selected, the backlight setting is determined by the selected LUT.
		If "Disabled" is selected, the "Backlight" setting in the "Picture" menu applies. Note: The backlight retains this value even if the LUT is changed.
LUT	LUT 1-5,	Select the Look Up Table (LUT)
	LUT User	The LUT determines the monitor's gamma
	lows s	curve. Using different LUTs for example, allows specific grayscale levels to be highlighted.
		When the "User" LUT is selected, the LUT can be individually set through the "Adjust LUT" OSD menu. See also "Adjust LUT" menu [> 51].
		A table with information regarding the available LUTs is displayed below the LUT selection.
		Note: Select a DICOM LUT to view radiographic images.

7.10 "Info" menu

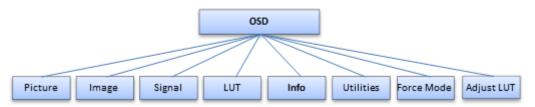

Function	Values	Description
Information	For display only	Displays the current monitor and operating data in the respective picture mode
		Serial number
		Temperature [°C]
		Operation time [h:m]
		Backlight operation time [h:m]
		BL on over 350 cd/m² [h]
		Source
		Mode
		Source signal
		2nd Source
		Mode
		2nd Source signal
		FW type
		FW version
		OSD version
		Config version
		SDK version

7.11 "Utilities" menu

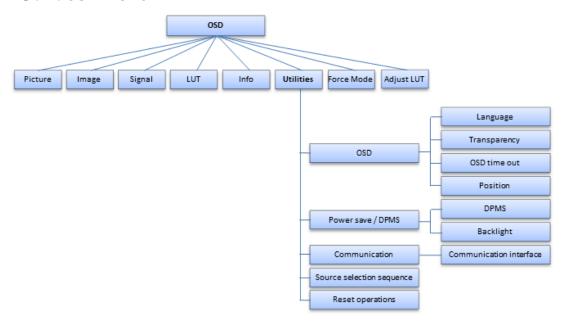


Function	Values	Description
OSD		
Language	German	Set language
	English	German or English can be selected.
	Default: English	
Transparency	Opaque	OSD background selection.
	Transparent	
	Half transparent	
	Default: Opaque	
OSD timeout	Off / 5 sec / 10 sec / 30 sec / 1 min	Closes the OSD menu after the specified time of inactivity.
	Default: Disabled	
Position	Top - left / Top - center / Top - right / Mid - left / Centered / Mid - right / Bottom - left / Bottom - center / Bottom - right	Set the location of the OSD menu on the screen area using predefined positions
	Default: Bottom - right	

Function	Values	Description
Power save / DPMS		
DPMS	Enabled / Disabled Default: Enabled	Enabling or disabling the DPMS (Display Power Management System)
	Dordan: Endord	Note: "Enabled" and "Disabled" modes of the DPMS standard are supported for analog signals. For digital signals, DVI-DMPM is supported at the DVI connection sockets
		If "DPMS" is activated and there is no input signal, the backlight is switched off or dimmed. This saves power, and increases the service life of the backlight.
Backlight	Dimmed Disabled	The "Backlight" can be set to "Dimmed" or "Off"
	Default: Disabled	When the DPMS mode is active, the backlight is either dimmed or switched off. The controller board and power supply unit remain active. For digital signals, 5 V is also applied to the DVI socket for external modules.
Communication		
Communications inter-	DDC/DVI 1 digital	Select the interface for communication
face	DDC/DVI 2 digital	Serial communication is always active.
	USB	Installers can also select between USB and
	Default: USB	DDC.
Source selection sequence		The source scan is interrupted as soon as the first source with an active video signal is found. This signal is then displayed.
Reset operations	Reset to firmware de- faults	All user settings are deleted and reset to the firmware defaults.

7.12 "Force Mode" menu

Force Mode is used to set the timing of video sources that cannot be displayed with the timings saved in the monitor.

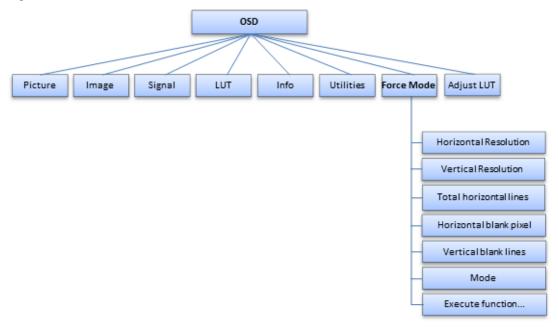


Function	Values	Description	
Horizontal Resolution		Set the horizontal resolution in the active resolution range.	
Vertical Resolution		Set the vertical resolution in the active resolution range.	
Total horizontal lines		If the resolution is known, this can be used as the approximation value. The following applies: "Total horizontal lines" > "Horizontal resolution".	
		A larger "Horizontal resolution" can result in the image being split vertically.	
Horizontal blank pixel		Set the horizontal blanking interval.	
Vertical blank pixel		Set the vertical blanking interval.	
Mode	Active / Start /		
	Inactive Default: Inactive	When "Active" is selected, the specified timing values are applied with "Execute function".	
		Force Mode is switched off when "Inactive" is selected.	
		When "Start" is selected, the start values for the timing are determined using "Execute function".	
		See also Determining the timing data	
Execute function		When "Start" is selected in Mode, the set video resolution is executed.	

7.13 "Adjust LUT" menu

Note

- The "Adjust LUT" menu is displayed only if the "User" setting is selected in the "LUT"
- After each parameter change, the LUT in the device is immediately recalculated. This is indicated in the OSD menu by the message: "Gamma table being recalculated". When the calculation is complete, the message disappears from the OSD.

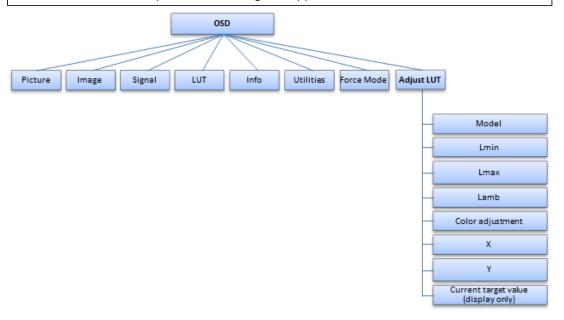

/!\CAUTION

Check the gray levels

After each LUT recalculation performed in the "Adjust LUT" menu, it is absolutely necessary to apply an appropriate test pattern and check that all gray levels are displayed correctly.

Function	Values	Description
Model	NativeCIE 1976	The gamma model, which serves as the basis for re- calculating the LUT, can be selected here.
	 CIE (DIN 6174) CRT DICOM Gamma 2.0, 2.2, 2.5, 2.7, 3.0, 3.3 Log.Lum.linear 	Note: When "Natural" is selected, only the value for Lmax can be changed. All other functions are hidden.
	Default: DICOM	
Lmin	0 - 250	Adapt the minimum black value between 0 and 5 cd/
	Default: 40 (0.8 cd/m ²)	m ²

Operation 7.13 "Adjust LUT" menu

Function	Values	Description
Lmax	0 - 250 Default: 80 (400 cd/m²)	Adapt the maximum white value between 200 and 600 cd/m ²
Lamb	0 - 250 Default: 25 (0.5 cd/m²)	Takes into account the ambient light between 0 and 5 cd/m² when recalculating the LUT.
Color adjustment	BluishIntrinsicTypicalDefault: Intrinsic	 Set the color coordinates to the following predefined values: Bluish: x = 0.272; y = 0.308. Intrinsic: Color coordinates remain unchanged from their intrinsic values. Typical: x = 0.300; y = 0.315.
х	-100 - +100 Default: 0 (x = 0.272)	X color coordinate correction When using the color setting based on the existing defaults, it is also possible to individually modify the color coordinates. The X color coordinate can be adapted using the "X" function.
Υ	-100 - +100 Default: 0 (y = 0.308)	Y color coordinate correction When using the color setting based on the existing defaults, it is also possible to individually modify the color coordinates. The Y color coordinate can be adapted using the "Y" function.
Current target value (display only)		

8 Cleaning and check settings

8.1 Cleaning

Recommended cleaning agents and disinfectants

NOTICE

Device maintenance, cleaning and disinfecting

- Make sure liquids do not seep into the device. Liquids that seep into the device may result in an electric shock or failure of the device.
- The screen is extremely sensitive to mechanical influences. Absolutely avoid scratches, shocks, etc. for this reason.
- Clean the screen when dirty using a microfiber cloth and, if necessary, a recommended cleaning agent. Clean the housing parts with a recommended cleaning agent.
- Use only tested disinfectants.
- If a cleaning agent is sprayed directly onto the screen surface, use a microfiber cloth to remove drops which run down before they reach the edge of the panel.
- Remove drops of liquid from the device immediately. Contact with liquids over a longer period can cause discoloration or allow calcium deposits to form on the surface.

Agent class	Tested cleaning agents and disinfectants:	Further examples
Aldehyde	Melsitt	Aldasan 2000
		Kohsolin
		Gigasept FF
		Cidex
Alcohol	Ethyl alcohol, 96%	Hospiset cloth
Peroxide compounds	Perform	Dismozon pur
Chlorine derivatives	Terralin	Quartamon Med
Disinfectants	Mikrozid liquid	TaskiDS5001 (Diverseylever Labs)
	Nocolyse	Morning Mist
		Surfanios Fraicheur Citron (Anios Labs)
Guanidine derivatives	Lysoformin	
Quaternary compounds	Incidur spray, undiluted	
Standard household washing-up liquid	denk mit	Fairy Ultra, Pril, Palmolive
Pyridine derivatives	Activ spray, undiluted	
Water	Tap water	
	Distilled water	

/!\CAUTION

• Information on cleaning or disinfection of other system components can be obtained from the respective instructions for use.

Prohibited cleaning agents and disinfectants

The following cleaning agents and disinfectants can bleach the paint after a longer period of application:

Agent class	Tested cleaning agents and disinfectants:	Further examples
Petroleum spirit	Petroleum spirit close to boiling	Petroleum ether

8.2 Check the settings

/ CAUTION

Checking the settings

- The settings may only be checked by service personnel.
- The settings must not be checked in the presence of patients.

The picture quality of the monitor changes due to aging of the LCD unit and the backlight.

- Check the monitor settings at regular intervals in accordance with the local guidelines.
- Correct the settings if necessary.

9 Troubleshooting

In normal operation, an image is displayed when the sources are set correctly. In the event of a fault, the following checks can help localize the cause, based on the screen display and the signal information in the OSD.

- 1. Check the device for the possible causes listed in the following.
- 2. Carry out the remedial measures before contacting service personnel.

9.1 Fault correction

Fault	Cause	Remedy
No picture appears on the monitor	Monitor is in standby	Check whether the monitor is in energy saving mode; for example, has Power Management been activated by the signal source and monitor?
		End standby by holding down "Standby" on the remote control.
	Power cable not connected	Check the power cable.
	or connected incorrectly	Is the power switch in the "On" position?
		Can the OSD menu be activated?
	No sync signal or video signal	Check whether source sweep has been switched on and the source chosen under source selection is connected.
		Check the video cable.
	Sync cable connected incorrectly	Exchange the sync cable connectors.
	The video source is not	Check the video source.
	sending a video signal, or is sending an unknown timing.	In the "Info" OSD menu, check whether the source was detected and that the source signal is in the displayable range.
	The video source is sending an unknown timing	Set a different timing for the video source.
	Blown fuse	Contact service.
Fuzzy picture, interference in vertical lines	Scanning frequency or phase incorrectly set.	Set the frequency and phase.
Picture is without contrast and has a strong	The video source is only transmitting a green signal	In the OSD "Signal" menu, switch from "Color" to "Monochrome".
green shading		Check the selected LUT and whether the USER settings can be changed.
		Has the device been running for at least 20 minutes without Power Management activated?
Other	Loose plug	Connect the loose plug and secure it.
	Defective cable	Replace the defective cable.

10 Technical specifications

Note

Applicability of technical specifications

All technical specifications are valid after a warm-up period of 30 minutes.

10.1 Monitor characteristics

Display type	Color, TFT, IPS
Active area	1073.8 mm x 604 mm
Screen diagonal	48.5" (123.2 cm)
Resolution	1920 x 1080 pixels
Refresh rate	60 Hz
Pixel arrangement	RGB vertical stripes
Pixel distance	0.55926 mm x 0.55926 mm
Contrast ratio	1300:1 typical, 900:1 minimal
Horizontal viewing angle	Typically 178° for contrast ratio > 10:1
Vertical viewing angle	Typically 178° for contrast ratio > 10:1
Backlight type	White LEDs
Screen brightness	Typically 700 cd/m ²
	Min. 560 cd/m ²

10.2 Power supply

Power connector	Non-heating appliance socket
Line voltage	AC 100 - 240 V AC
Miniature fuse	2 T6.3 AH: UL 248-14, 250 V (5 x 20 mm)
Maximum current consumption	1.5 - 0.7 A
Power consumption	85 W (typical), 144 W (max)
Power consumption in Standby Mode	19 W

10.3 Electronics

Multi-standard technology	Video modes with resolutions less than 1920 x 1080 can be expanded to the TFT resolution, and thus utilize the full display area (like multi-sync CRTs)	
	For images with a resolution other than 1920 x 1080, select "One to one" for "Image size / zoom" in the "Image" menu. Otherwise, you will have to deal with scaling artifacts.	
	Caution: If the timing is frame-buffered or frame-synchro- nized, the picture information may be lost; the gray scales - the color hues for color images - are also reduced and may be invisible.	
Timing recognition	The monitor detects various signals at the inputs and automatically adjusts itself to timings such as the following:	
	VGA: 60-85 Hz	
	• SVGA: 56-85 Hz	
	• XGA: 60-85 Hz	
	• SXGA: 60-85 Hz	
	• FHD: 60 Hz	
	UXGA: 60 Hz with vertical section	
	PAL: RGB, YUV	
	NTSC: RGB, YUV	
	• SD-SDI: 487i, 576i	
	HD-SDI: 720p, 1080i, 1080p30	
	• 3G-SDI: 720p, 1080i, 1080p30, 1080p60	
Communication interface	RS-232: Via RJ 11 socket	
	USB 2.0: 2 downstream (type A), 1 upstream (type B)	
	DDC: Via DVI-I	

10.4 Mechanical design

Housing components	Metal	
Housing color	• White (NCS S 1000-N)	
	Black (RAL 7021)	
Protective glass	Optically non-reflective protective glass on both sides	
Ventilation openings	In rear panel	
Degree of protection according to	Front: IP65	
EN 60529	Back: IP20	
Connection panel	On rear panel, under cover	
Weight in kg	31 kg +/- 10%	
Dimensions (W x H x D) in mm	1103 x 634 x 80.9	

10.5 Climatic conditions

In operation	
Temperature range	+5 °C - +40 °C ambient temperature
Temperature gradient	Max. 5 °C/h, no condensation
Air pressure	700 - 1060 hPa (3000 m360 m)

For transport and storage (packed)	
Temperature range	-20 °C - +60 °C ambient temperature
Temperature gradient	Max. 5 °C/h, no condensation
Humidity	10 - 90 %, non condensing, at 25 °C
Air pressure	200 - 1060 hPa (12000 m360 m)


10.6 Safety regulations

Safety standards	• IEC 60601-1			
	• IEC 60950-1			
	• CAN/CSA - C 22.2 60601-1			
	GB 4943.1 (non-tropical, altitude < 2000 m)			
Protection class	Protection class I			
Degree of protection according to	Front: IP65			
EN 60529	Back: IP20			
Medical device classification (EU)	Class 1			


11 Dimension drawings

All dimensions in mm

11.1 View from front and side

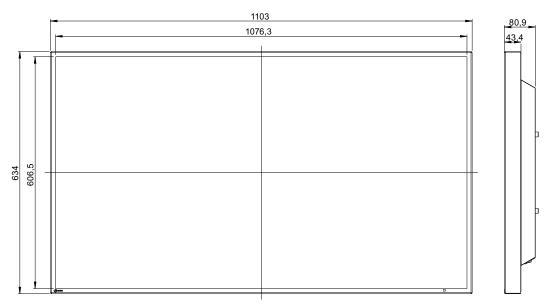

11.2 Rear view

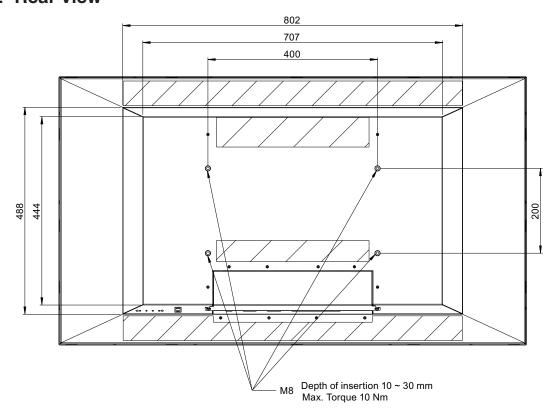

12 Spare parts / accessories

12.1 Accessories

Accessories	Order number	Description
FST4700 Black	6GF69888BA03	Stand, black
FST4700 White	6GF69888BA04	Stand, white

13 Appendix

13.1 Information on electromagnetic compatibility (EMC)

EIZO monitors were designed for the display of images and normal monitor operation.

/ WARNING

Special EMC provisions are required for use of the CuratOR LX491W. Installation, assembly, and use must take place in compliance with the following instructions.

- Only use the cables included in the scope of delivery or recommended by the manufacturer. The use of other cables can result in increased electromagnetic radiation and reduced electromagnetic interference immunity of the device, as well as improper use.
 Cable length: max. 3 m
- The monitor should not be placed on other devices or positioned in their immediate vicinity. If devices have to be operated on or in the immediate vicinity of one another, the monitor or system must be monitored to ensure proper operation for the defined configuration.
- When using a portable RF communications device, maintain a distance of at least 30
 cm from all parts of the monitor, including cables. Otherwise, problem-free function of
 the device cannot be guaranteed.
- Persons connecting additional devices to the signal input or output when configuring a medical system are responsible for ensuring compliance with standard IEC/EN 60601-1-2.

Electromagnetic radiation

The CuratOR LX491W is intended for use in the electromagnetic environments noted below.

Customers and users of the CuratOR LX491W have to ensure that the device is used in such an en-

Customers and users of the CuratOR LX491W have to ensure that the device is used in such an environment.

Radiation test

Conformity

Information regarding the electromagnetic environment

Radiation test	Conformity	information regarding the electromagnetic environment
RF radiation CISPR11/EN 55011	Group 1	The CuratOR LX491W generates RF for its internal function only. For this reason, the RF radiation is very low and is therefore unlikely that it will result in interference in electronic devices in the immediate vicinity.
RF radiation	Class B	The CuratOR LX491W is approved for use in a number of en-
CISPR11/EN 55011 GB9254		vironments. This includes residential areas and areas connected directly to the public low-voltage grid, such as private homes.
Harmonic currents	Class D	
IEC/EN 61000-3-2 GB17625.1		
Voltage fluctuations / flicker	fulfilled	
IEC/EN 61000-3-3		

Electromagnetic interference immunity

The CuratOR LX491W was tested with the following compliance levels in accordance with the test requirements for professional healthcare facilities, as established in IEC/EN 60610-1-2.

Customers and users of the CuratOR LX491W have to ensure that the monitor is used in such an environment.

Interference im- munity test	Measurement level	Compliance level	Information regarding the electromagnetic environment
Electrostatic discharge (ESD) IEC/EN 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	It is recommended to use the device on wood, concrete, or ceramic floors. If the floor is made of synthetic material, the relative humidity should be at least 30%.
Fast transient electric distur- bances (bursts) IEC/EN 61000-4-4	±2 kV power lines ±1 kV input / out- put lines	±2 kV power lines ±1 kV input / out- put lines	The power supply quality has to correspond to that of typical industrial environments or hospitals.
Surge voltage IEC/EN 61000-4-5	±1 kV line against line ±2 kV line against ground	±1 kV line against line ±2 kV line against ground	The power supply quality has to correspond to that of typical industrial environments or hospitals.
Voltage dips, brief interrup- tions, and fluctu- ations of power supply lines IEC/EN 61000-4-11	0 % V_T for 0.5 periods and 1 period 70 % V_T for 25 / 30 periods at 50 / 60 Hz 0 % V_T for 250 / 300 periods at 50 / 60Hz	0 % V_T for 0.5 periods and 1 period 70 % V_T for 25 periods at 50 Hz 0 % V_T for 250 periods at 50 Hz	The power supply quality has to correspond to that of typical industrial environments or hospitals. If the monitor has to continue operation even if the power supply is interrupted, it is recommended to connect the device to an uninterruptible power supply or battery.
Magnetic fields with energy technology fre- quencies IEC/EN 61000-4-8	30 A/m (50 / 60 Hz)	30 A/m (50 Hz)	The magnetic fields with energy technology frequencies must be in an area that is representative of a typical location in a typical industrial environment or hospitals. The device should be used at least 15 cm away from the source of magnetic fields with energy technology frequencies.

Note: V_T is the alternating current voltage before application of the measurement level.

Electromagnetic interference immunity

The CuratOR LX491W was tested with the following compliance levels in accordance with the test requirements for professional healthcare facilities, as established in IEC/EN 60601-1-2. Customers and users of the monitor have to ensure that the monitor is used in such an environment.

Interference immunity test	Measure- ment level	Compliance level	Information regarding the electromagnetic environment
Line-based dis- turbances caused by RF	3 V _{rms} 150 kHz to 80 MHz	3 V _{rms}	Portable and mobile RF communications devices may only be operated in the vicinity of the monitor and its components (including cables) when in com-
fields IEC/EN 61000-4-6	6 V _{rms} ISM bands between	6 V _{rms}	pliance with the recommended minimum distance. It is determined using the formula for calculating the frequency of the transmitter.
	150 kHz and		Recommended minimum distance
	80 MHz		$d = 3.5/3 \ \sqrt{P} = 1.2 \ \sqrt{P}$, 150 kHz to 80 MHz
Electromagnetic	•		d = 2 √P, ISM bands between 150 kHz and 80 MHz
RF fields	80 MHz to 2.7 GHz	*=	$d = 3.5/3 \ \sqrt{P} = 1.2 \ \sqrt{P}$, 80 MHz to 800 MHz
IEC/EN 61000-4-3			$d = 7/3 \ \sqrt{P} = 2.3 \ \sqrt{P}$, 800 MHz to 2.7 GHz
61000-4-3			In this case, "P" stands for the measured maximum nominal output power in watts (W) of the transmitter recommended by the transmitter manufacturer, and "d" for the recommended minimum distance in meters (m).
			The field strengths of fixed transmitters according to electromagnetic location measurement ^{a)} have to be less than the compliance level in each individual frequency range.
			Interference can occur when used in the vicinity of devices identified with the following symbol.

Note: The higher frequency range applies at 80 MHz and 800 MHz.

Note: Guidelines with respect to line-based interference due to RF fields or electromagnetic RF fields may not apply in all situations. Absorption and reflection by structures, objects, and people impact the propagation of electromagnetic waves. .

^{a)} The field strengths of fixed transmitters, for example the base station for cordless and mobile telephones, radio, land mobile radio, ham radio, and television cannot be determined precisely in advance. To evaluate the electromagnetic environment using fixed transmitters, an electromagnetic location measurement should be included. If the measured field strength in the environment where the device is used exceeds the applicable RF compliance level, observe the monitor to ensure its proper operation. If improper operation is observed, in some circumstances additional measures may be necessary, such as reorienting or repositioning the device.

13.1 Information on electromagnetic compatibility (EMC)

Recommended minimum distance between portable or mobile RF communications devices and the CuratOR LX491W

The CuratOR LX491W is intended for use in an electromagnetic environment in which interference due to electromagnetic radiation is controlled. For other portable and mobile RF communication devices (transmitters), the recommended minimum distance between the portable and mobile RF communication devices (transmitters) and the monitor applies as listed below. This is based on the maximum output power of the communication device.

Maximum nominal output power of	Recommended minimum distance according to the frequency of the transmitter (m)					
the transmitter (W)	150 kHz to 80 MHz d = 1.2 √P	80 MHz to 800 MHz d = 1.2 √P	800 MHz to 2.7 GHz d = 2.3 √P			
0.01	0.12	0.12	0.23			
0.1	0.38	0.38	0.73			
1	1.2	1.2	2.3			
10	3.8	3.8	7.3			
100	12	12	23			

For transmitters whose maximum nominal output power is not shown above, the recommended minimum distance "d" in meters (m) can be determined using the formula for calculating the frequency of the transmitter. "P" here stands for the transmitter's maximum measured nominal output power in watts (W), as recommended by the transmitter's manufacturer.

Note: For 80 MHz and 800 MHz, the recommended minimum distance for the higher frequency range applies.

Note: This information may not be applicable in all situations. Absorption and reflection by structures, objects, and people impact the propagation of electromagnetic waves.

Recommended minimum distance between portable or mobile RF communications devices and the CuratOR LX491W

The CuratOR LX491W is intended for use in an electromagnetic environment in which interference due to electromagnetic radiation is controlled. The customer or user of the monitor can help prevent electromagnetic interference by maintaining the recommended minimum distance between portable and mobile RF communications devices (transmitters) and the monitor.

The interference immunity regarding adjacent fields has been confirmed for the following wireless RF communications devices:

Test frequency (MHz)	Band- width ^{a)} (MHz)	Service ^{a)}	Modula- tion ^{b)}	Maximum power (W)	Minimum distance (m)	Measure- ment level (V/m)	Compli- ance level (V/m)	
385	380 - 390	TETRA 400	Pulse mod- ulation ^{b)} 18 Hz	1.8	0.3	27	27	
450	430 - 470	GMRS 460 FRS 460	FM ±5 kHz deviation 1 kHz sine	2	0.3	28	28	
710	704 - 787	LTE band 13,	Pulse mod-	0.2	0.3	9	9	
745		17	ulation ^{b)} 217 Hz					
780								
810	800 - 960	GSM 800/900	Pulse mod-	2	0.3	28	28	
870		TETRA 800 iDEN 820	ulation ^{b)} 18 Hz					
930		CDMA 850 LTE band 5						
1720	1700 -	CDMA 1900 ulation ^b	Pulse mod-	ation ^{b)}	0.3	28	28	
1845	1990		217 Hz					
1970		DECT LTE band 1, 3, 4, 25 UMTS						
2450	2400 - 2570	Bluetooth WLAN 802.11 b/g/n RFID 2450 LTE band 7	Pulse mod- ulation ^{b)} 217 Hz	2	0.3	28	28	
5240	5100 -	WLAN 802.11	Pulse mod-	0.2	0.3	9	9	
5500	5800	a/n		ulation ^{b)} 217 Hz				
5785								

^{a)} For some radio services, only the frequencies for the radio contact from the mobile communications device to the base station ("Uplink") is included in the table.

b) The carrier is modulated by a square wave with 50 % duty cycle.

13.2 Markings and symbols

Marking/symbol	Description
\triangle	Symbol for "Caution, observe accompanying documents".
4	Symbol for "Dangerous voltage".
CE	CE marking (EU conformity mark).
EU Medical Device	Medical device in accordance with the European medical device regulation.
Electrical Safety WET o us E113208	MET marking, in accordance with U.S. and Canadian national regulations.
(1)	CCC marking, in accordance with for Chinese national regulations (pending).
	RCM marking for conformity with Australian EMC standards (pending).
F©	U.S. FCC marking for communication devices.
	Symbol for manufacturer of medical devices, supplemented by the date of manufacture.
X	WEEE marking: Product must be disposed of separately; materials may be recycled.
©	Marking according to ACPEIP (China-RoHS).
IP20	Symbol for International Protection Marking according to IEC 60529.
	Symbol for protective ground (ground)
	"On" symbol (voltage)
	"Off" symbol (voltage)
	Input for service calls
•~	Symbol for USB
⊕	Symbol for input.
\rightarrow	Symbol for output.
===	Symbol for direct current.
	Symbol "Comply with the instructions for use".

13.3 FCC Declaration of Conformity

For U.S.A., Canada, etc. (rated 100-120 Vac) Only

FCC Declaration of Conformity

We, the Responsible Party

EIZO Inc.

5710 Warland Drive, Cypress, CA 90630

Phone: +1 (562) 4 31 50 11

declare that the productTrade name: EIZO

Model: CuratOR LX491W

is in conformity with Part 15 of the FCC Rules. Operation of this product is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures.

- Reorient or relocate the receiving antenna.
- · Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note

Use the attached specified cable below or EIZO signal cable with this monitor so as to keep interference within the limits of a Class B digital device.

- AC Cord
- Shielded Signal Cable (enclosed)

Canadian Notice

This Class B digital apparatus complies with Canadian ICES-003.

Cet appareil numérique de la classe B est comforme à la norme NMB-003 du Canada.

13.4 China RoHS (Restriction of Hazardous Substances)

13.4 China RoHS (Restriction of Hazardous Substances)

液晶显示器 LCD Monitor

型号 Model: CuratOR LX491W

根据SJ/T11364-2014《电子电气产品有害物质限制使用标识要求》特提供如下有关污染控制方面的信息。

The following product pollution control information is provided according to SJ/T11364-2014 Marking for the restriction of the use of hazardous substances in electrical and electronic product.

电子信息产品污染控制标志说明 Explanation of Pollution Control Label



该标志表明本产品不含有超过中国标准GB/T26572-2011《电子电气产品中限用物质的限量要求》中限量的有毒有害物质。该标志还表示本产品废弃后可以回收利用,不应随意丢弃。

This symbol indicates the product does not contain hazardous materials in excess of the limits established by the Chinese standard GB/T26572-2011 Requirements of concentration limits for certain restricted substances in electrical and electronic products. The symbol also signifies that the product can be recycled after being discarded, and should not be casually discarded.

有毒有害物质或元素的名称及含量 Name and Concentration of Hazardous Substances

部件名称 Component Name	有毒有害物质或元素 Hazardous substances' name					
	铅 (Pb)	汞 (Hg)	镉 (Cd)	六价铬 (Cr(VI))	多溴联苯 (PBB)	多溴二苯醚 (PBDE)
液晶纯平屏幕 LCD Flat Screen	0	0	0	0	0	0
控制板 Controller Board	0	0	0	0	0	0
电源 Power Supply	0	0	0	0	0	0
其他 电路板 Other Circuit Boards	0	0	0	0	0	0
其他(电缆等) Others (cables, etc.)	0	0	0	0	0	0
机架、底盘 Housing, Chassis	0	0	0	0	0	0
附件(信号电缆、输电线等) Acessories (signal cable, power line, etc.)	0	0	0	0	0	0

部件名称 Component Name

有毒有害物质或元素 Hazardous substances' name

本表格依据SJ/T 11364 的规定编制。

- O: 表示该有害物质在该部件所有均质材料中的含量均在GB/T 26572 标准规定的限量要求以下
- X: 表示该有害物质至少在该部件的某一均质材料中的含量超出GB/T 26572 标准规定的限量要求
- 此表所列数据为发布时所能获得的最佳信息.
- 由于缺少经济上或技术上合理可行的替代物质或方案,此医疗设备运用以上一些有害物质来实现设备的预期临床功能,或给人员或环境提供更好的保护效果。

This list is based on SJ/T 11364.

- O: Indicates that this toxic or hazardous substance contained in all of the homogeneous materials for this part is below the limit requirement in GB/T 26572.
- X: Indicates that this toxic or hazardous substance contained in at least one of the homogeneous materials used for this part is above the limit requirement in GB/T 26572.
- Data listed in the table represents best information available at the time of publication.
- Applications of hazardous substances in this medical device are required to achieve its intended clinical uses, and/or to provide better protection to human beings and/or to environment, due to lack of reasonably (economically or technically) available substitutes.

产品中有毒有害物质或元素的名称及含量 Table of hazardous substances' name and concentration.

13.5 Declaration of compliance with India RoHS

We, EIZO Corporation, hereby declare and guarantee that this product has been designed and manufactured in compliance with the E-Waste management rule 2016 which prohibit the inclusion of the following substances except for the exemptions listed in schedule II.

- Lead, Mercury, Hexavalent Chromium, Polybrominated Biphenyls or Polybrominated
 Diphenyl Ethers exceeding a concentration of 0.1% by weight in homogeneous materials
- Cadmium exceeding a concentration of 0.01% by weight in homogeneous materials

For information on proper disposal and recycling of the product, please refer to the following website.

eizo.co.in/e-waste.php

13.6 Environmental protection

Comply with all local requirements and laws pertaining to the disposal of devices.

The device is in compliance with directive 2011/65/EU for limiting the use of specific hazardous materials in electric and electronic devices.

13.7 Additional devices

Connected devices must meet the relevant safety standards.

13.8 Warranty

Opening of the housing, or electrical or mechanical changes on or in the device, result in cancellation of the warranty. For warranty details, please contact the sales partner from whom you purchased the product. These warranty conditions are neither extended nor limited by the contents of this instruction manual.

13.9 Repairs

13.9 Repairs

Please contact the sales partner from whom you purchased the product.

13.10 Contact

Support during installation and for technical questions

www.eizo-or.com

13.11 Trademarks

The EIZO Logo is a registered trademark of EIZO Corporation in Japan and other countries.

EIZO is a registered trademark of EIZO Corporation in Japan and other countries.

RadiForce is a registered trademark of EIZO Corporation in Japan and other countries.

CuratOR is a registered trademark of the EIZO Corporation.

RadiCS is a registered trademark of EIZO Corporation in Japan and other countries.

RadiNET is a registered trademark of EIZO Corporation in Japan and other countries.

ScreenManager is a registered trademark of EIZO Corporation in Japan and other countries.

Windows is a registered trademark of Microsoft Corporation in the United States and other countries.

Apple is a registered trademark of Apple Inc.

Macintosh is a registered trademark of Apple Inc.

Mac is a registered trademark of Apple Inc.

VESA is a registered trademark of the Video Electronics Standards Association in the United States and other countries.

HDMI is a trademark or registered trademark of HDMI Licensing, LLC in the United States and other countries.

DICOM is a registered trademark of the National Electrical Manufacturers Association for the publication of standards in the area of digital communication of medical information.

TORX is a registered trademark of Acument Intellectual Properties, LLC.

All other company and product names are trademarks or registered trademarks of their respective owners.

Index

A		0	
Assembly	16	OSD menu Adjust LUT Force Mode Image quality	51 50 37
Cleaning staff	6	Image quality Image setting	39
Connecting	21	Info	47
Connection panel	20	LUT	46
Contact	70	Overview	37
		Signal	42
В		Tools	48
D		Output	20
Disposal	69	Overheating	14
E		Р	
Energy saving	23, 24	Pixel defects	24
6,7 6	•	Power connector	22
F		Protective screen	13
Fault	55		
Force Mode	26, 50	S	
		Service personnel	6
•		Shielding measures	17
G		Signal input	20, 21
General safety instructions	7	Signal output	20
Grounding	17		
1		T	
1		Technical specifications	56
Image Sticking	23	Timing	
Input	20	Basic information	27
Installation site	14	determining	29
		Troubleshooting	55
<u>L</u>			
Line power	22	U	
LUT	46	User	6
Adjust	51		
М		W	
Monitor information	47	Warranty	69
wonto incinduon	+ /		



EIZO GmbH

Carl-Benz-Straße 3 76761 Rülzheim Germany